


# Dual Applicator Quick Start Guide

Feature only available on the A2500 PRO device

Pulsed Harmonix now offers the ability to operate two accessory applicators from a single True Pulse **A2500 PRO** PEMF wellness device!

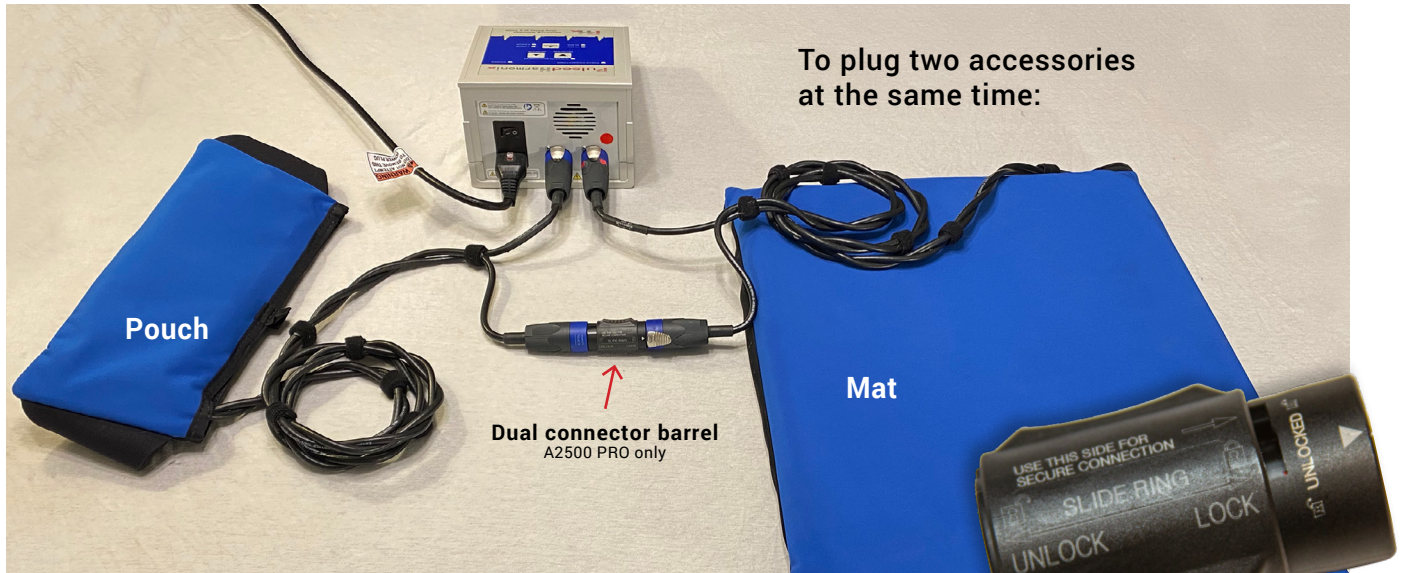

To plug two accessories at the same time:



**1.** Insert and twist the red dot connector of one applicator into the red dot socket on the Pulsar. Insert and twist the non-red dot connector of the other applicator into non-red dot socket on the Pulsar.



**2.** Insert either remaining connector into the Barrel and twist to secure. Insert the last connector into the Barrel and twist to secure.



**3.** Check that all four connectors are locked in by their slide tab.



**4.** The circuit is now complete.